

BloodSTOP®
CONTROLS BLEEDING FAST

*Power to control bleeding
Power to heal*

BloodSTOP® **Wound Care Dressing** *Water Soluble Hemostat*

“With BloodSTOP®, hemostasis
is achieved in seconds.”

- David Ngo, RN
King Drew Medical Center
Orange County, CA

Athletes

Schools

Home Healthcare

General First Aid

Paramedic Emergency

LifeScience PLUS
REVOLUTIONARY WOUND CARE

Rev 2018/03

BloodSTOP®

STOPS BLEEDING FAST




BloodSTOP® is a next generation hemostatic dressing made from 100% purified, water soluble, etherified sodium carboxymethyl cellulose. It is a woven matrix of natural plant cellulose, used by medical professionals and hospitals. When applied to a wound, BloodSTOP® quickly absorbs blood and other body fluids, transforms into a translucent gel to seal the wound and creates a moist environment for optimal natural wound healing. BloodSTOP® exerts its pro-coagulant activity by activating the intrinsic coagulation pathway, accelerating clot formation.

BloodSTOP® has EU class IIb CE certification for surface wound use, including with anticoagulant medications. In the US, BloodSTOP® has FDA 510K market clearance for surface wound hemostasis.

Features


- Acts fast
- Conforms to all wound surfaces
- 100% natural component
- Active bleeding control over traditional gauzes
- Painless upon application and removal
- Water soluble for easy removal without disruption to wound surface
- Adheres to the wound without the discomfort of adhesives
- Stays in wound while excess material and contaminants are easily washed off
- Easy to cut, pack or layer to fit any wound
- Biodegradable and biocompatible
- Cost effective

Instructions:

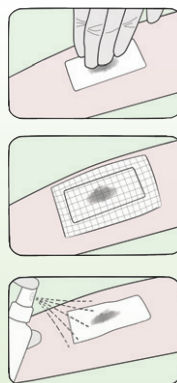

1. Apply BloodSTOP® and apply continuous pressure
2. Cover with second dressing
3. Remove with saline or water

For Nosebleeds:

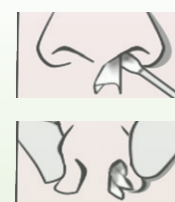

1. Insert
2. Compress
3. Rinse to remove



Product Description	Size	Quantity	Product Code
BloodSTOP® Hemostat	1"x1" (2.5x2.5cm)	20/Box	BS-09
	2"x2" (5x5cm)	20/Box	BS-10
	2"x4" (5x10cm)	20/Box	BS-11
	4"x4" (10x10cm)	10/Box	BS-12
	4"x8" (10x20cm)	10/Box	BS-13
	mini-pads 2x2cm	80/Box 4/Pk	BS-MP-18
	2"x2" (5x5cm) - 2 1"x1" (2.5x2.5cm) - 8	10/Box	BS-OTC-28